

WRIST


CLINICAL INDICATIONS: FX, Pain, swelling, Fibrocartilage tears , carpal tunnel										
parameter	SERIES 1	SERIES 2	SERIES 3	SERIES 4	SERIES 5	SERIES 6	SERIES 7	SERIES 8	SERIES 9	SERIES 10
description	3-PLN LOC	T1 AX	T2 AX FS	T1 SAG	GRE COR	STIR COR	T1 COR	3D T1 COR	T1 COR OR SAG FS C+	T1 AX FS C+
entry	Prone	Prone	Prone	Prone	Prone	Prone	Prone	Prone	Prone	Prone
position	head first	head first	head first	head first	head first	head first	head first	head first	head first	head first
coil	HRWRIST	HRWRIST	HRWRIST	HRWRIST	HRWRIST	HRWRIST	HRWRIST	HRWRIST	HRWRIST	HRWRIST
plane	3-PLN	AXIAL	AXIAL	SAGITTAL	CORONAL	CORONAL	CORONAL	CORONAL	COR OR SAG	AXIAL
pulse seq.	FGRE	FSE-XL	frFSE-XL	FSE-XL	GRE	STIR	FSE-XL	fSPGR	fSPGR	fSPGR
im options		EDR,TRF,VB	EDR,TRF,VB,Z512	VB,TRF		VB,TRF,Z512	VB,TRF	Z2,Z512	EDR,Z512	EDR,Z512
PSD Name										
Te		13	100	15	15	50	Min-Full	Min	OUT-OF-PHASE	OUT-OF-PHASE
TR		550	2850	50	500	3750	400	150	150	
TI						150				
flip angle					20			10	90	90
ETL		4	19	4	8	2				
RBW		16	20	15	11	16	12	11	31	31
SATS			FAT					FAT(special)	FAT(special)	FAT(special)
FOV	25	10	10	10	10	10	10	10	10	10
slice thick	5	3.5	3.5	3.5	3	3	3	1.8	3	3.5
spacing	5	0.5	0.5	0.3	0.3	0.3	0.3	0.3	0.3	0.5
F Phase	256	352	256	352	320	256	352	256	256	256
P Phase	128	162	192	192	192	192	224	192	192	192
NEX	1	3	4	3	2	3	3	3	3	3
Phase FOV	1	0.8	0.75	0.75	0.75	0.75	0.75	0.75	1	1
FREQ DIR		R/L	R/L	S/I	S/I	S/I	S/I	S/I	S/I	S/I
COVERAGE	WRIST	WRIST	WRIST	WRIST	WRIST	WRIST	WRIST	WRIST	WRIST	WRIST
contrast									1mmol/kg	1mmol/kg
comments:	IF CONTRAST INDICATED: RUN SERIES 9 AND 10									

WRIST


AXIAL



SAGITTAL

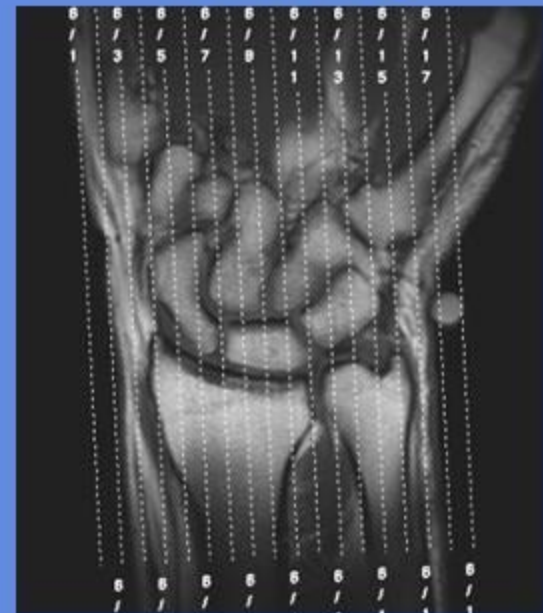

3D CORONAL

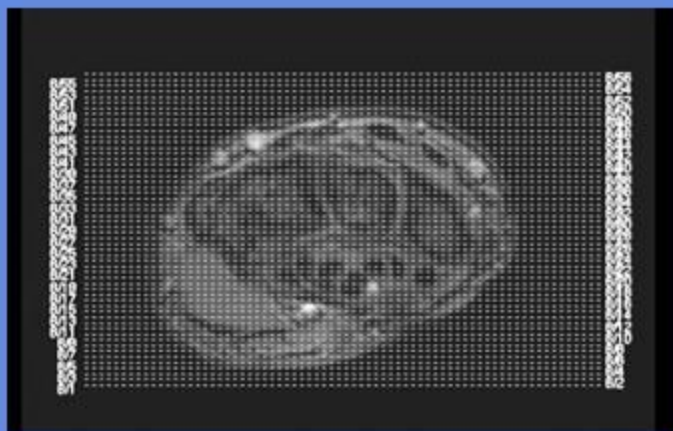

CORONAL

